

Annexin V-FITC/7-AAD Apoptosis Detection Kit

A1456530

Storage 2-8°C. Protect from light. Do not freeze.

Shipping Shipped with ice packs. Please store under the specified storage conditions immediately upon receipt.

Introduction

Annexins are a family of calcium-dependent phospholipid-binding proteins that preferentially bind phosphatidylserine (PS). Under normal physiologic conditions, PS is predominantly located in the inner leaflet of the plasma membrane. Upon initiation of apoptosis, PS loses its asymmetric distribution across the phospholipid bilayer and is translocated to the extracellular membrane leaflet marking cells as targets of phagocytosis. Once on the outer surface of the membrane, PS can be detected by fluorescently labeled Annexin V in a calcium-dependent manner.

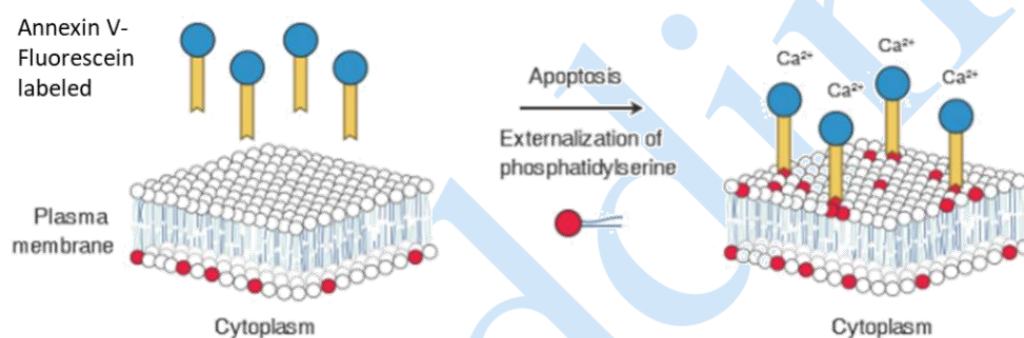


Fig. 1 Principal of Annexin-V Apoptosis Assay

In early-stage apoptosis, the plasma membrane excludes viability dyes such as propidium iodide (PI), 7-AAD. These cells will stain with Annexin V but not a viability dye, thus distinguishing cells in early apoptosis. However, in late stage apoptosis, the cell membrane loses integrity thereby allowing Annexin V to also access PS in the interior of the cell. A viability dye can be used to resolve these late-stage apoptotic and necrotic cells (Annexin V, viability dye-positive) from the early-stage apoptotic cells (Annexin V positive, viability dye-negative). This kit is suitable for the identification and enumeration of dead cells, such as apoptotic or necrotic cells, by flow cytometry.

Kit Contents

A1456530	Components	20T	50T	100T	Storage	Quantity Per Test
A1456530A	10×Annexin V Binding Buffer	5 mL	10 mL	20 mL	2-8°C.	200 µL per 0.5-1.0×10 ⁵ cells.
A1456530B	Annexin V-FITC	40 µL	100 µL	200 µL	2-8°C. Store in the dark.	2 µL per 0.5-1.0×10 ⁵ cells.
A1456530C	7-AAD Staining Solution	100 µL	250 µL	500 µL	2-8°C. Store in the dark.	5 µL per 0.5-1.0×10 ⁵ cells.

Note: The recommended number of cells for use is 0.5-1.0×10⁵ cells/test.

Instruction for use

1. Dilute 10× Binding Buffer to 1× using distilled water (1 mL 10× Binding Buffer + 9 mL ddH₂O).
2. Wash cells twice with cold PBS and then resuspend the desired amount of cells in Annexin V Binding Buffer at a concentration of 0.5-1.0×10⁶ cells /mL.
3. Add 2 µL of FITC Annexin V and 5 µL 7-AAD to 100 µL of the cell suspension.
4. Add 100 µL of 1× Binding Buffer to each assay. Gently vortex the cells and incubate for 10 min at RT (25°C) in the dark.
5. Analyze by flow cytometry within 1 hr.

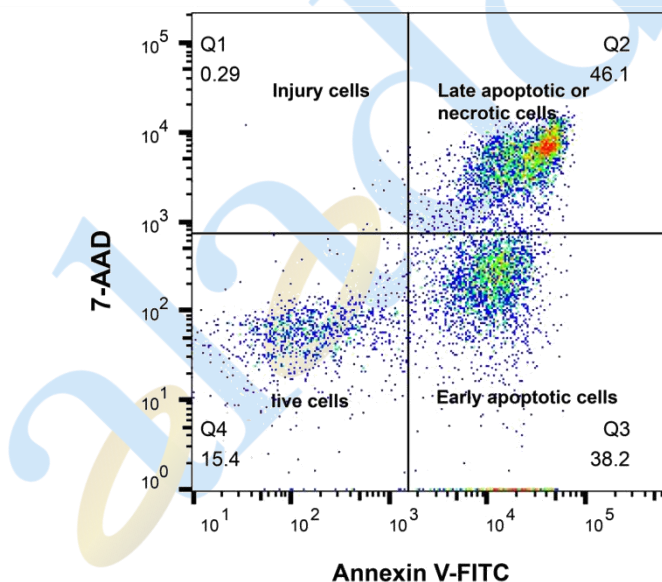


Fig. 2 Flow cytometry analysis of Jurkat cells treated with camptothecin using the Annexin V-FITC/7-AAD Apoptosis Detection Kit (A1456530). Q1: Injury cells, Loss of membrane integrity (Annexin V-) allows 7-AAD to enter the nucleus and bind DNA (7-AAD+). Q2: Late apoptotic or necrotic cells, Phosphatidylserine externalization (Annexin V+) with loss of

Web: <https://www.aladdinsci.com>

membrane integrity (7-AAD+). Q3: Live cells, Intact membrane without phosphatidylserine exposure (Annexin V-/7-AAD-). Q4: Early apoptotic cells, Phosphatidylserine externalization (Annexin V+) with intact membrane (7-AAD-).

Matters needing attention

1. Please try to avoid light when using to slow down the quenching of fluorescence.
2. 7-AAD solution is toxigenic and mutagenic; handle with care.
3. Due to the calcium dependence of the Annexin V:PS interaction, it is critical to avoid buffers containing EDTA or other calcium chelators during Annexin V experiments.

